ARDESTO™

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HEALTHY TEETH AND GUMS! EASY CLEANING OF ORTHODONTIC CONSTRUCTIONS!


You can effectively and easily take care of your gums, teeth, and beauty of your smile, using Ardesto irrigator on a daily basis! You are on your way to better gum health and a healthier smile! With the combination of water pressure and pulsation, the irrigator removes harmful bacteria and food debris deep between teeth and below gumline braces and other orthodontic constructions where brushing and even dental floss won't reach.

Traditional Clean

Traditional brushing can't easily reach hidden food debris and harmful bacteria



Flossing has limited access below the gumline and may damage gums

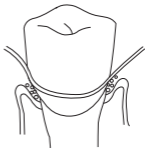Irrigator Clean

Removes trapped food debris and harmful bacteria with high frequency pulse water floss.

If you use irrigator on a daily basis, it gives you fresh breath, makes gums and teeth healthier, and helps to clean orthodontic constructions! Consult with your dentist on the irrigator modes usage.

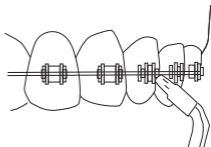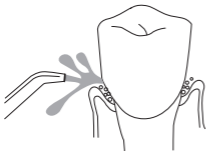CONTENT

Important Safeguards.....	21
Equipment.....	24
Product Structure.....	25
Product Description.....	26
Recommended technique.....	27
cleaning and Maintenance.....	32
Common Malfunction and Treatment.....	33
Correct disposal of this product.....	34
Warranty conditions.....	38

IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following.

READ ALL INSTRUCTIONS BEFORE USING

DANGER! To reduce the risk of electrocution:

- ◆ Always unplug this product immediately after using.
- ◆ Do not use while bathing.
- ◆ Do not place or store product where it can fall or be pulled into a tub or sink.
- ◆ Do not place in or drop into water or other liquid.
- ◆ Do not reach for a product that has fallen into water. Unplug immediately.

WARNING! To reduce the risk of

burns, electrocution, fire, or injury to persons.

- ◆ A product should never be left unattended when plugged in.
- ◆ This device can be used by children over 8 years of age and persons with physical or sensory or mental disabilities or with insufficient experience and knowledge, if they are under constant supervision or instructed in the safe use of the device and they understand the possible dangers. Cleaning and maintenance by the consumer should not be carried out by children without supervision.
- ◆ Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- ◆ Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into

water. Return the product to a service center for examination and repair.

- ◆ If the power cord is damaged, it must be replaced by the manufacturer, service center or similarly qualified personnel in order to avoid a hazard.
- ◆ Keep the cord away from heated surfaces.
- ◆ Never block the air openings of the product or place it on a soft surface, such as a bed or couch, where the air openings may be blocked. Keep the air openings free of lint, hair, and the like.
- ◆ Never use while sleeping or drowsy.
- ◆ Never drop or insert any object into any opening or hose.
- ◆ Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- ◆ Connect this product to a properly grounded outlet only.
See Grounding Instructions.


- ◆ Unplug this product before filling. Fill (reservoir) with water only unless otherwise specified by manufacturer. Do not overfill (or specify filling instructions).
- ◆ It's normal phenomenon that caused gum bleeding when you use the water flosser in the first time, which depended on different gum bearing ability.

SAVE THESE INSTRUCTIONS

EQUIPMENT

- Main Unit (Including Reservoir) x 1
- Instruction Manual (Warranty card included) x 1
- Ordinary High-Pressure Tip x 3
- Periodontal Tip x 1
- Orthodontic Tip x 1
- Dental Plaque Tip x 1
- Tongue Cleaner x 1
- Hydraulic rotating toothbrush x 1

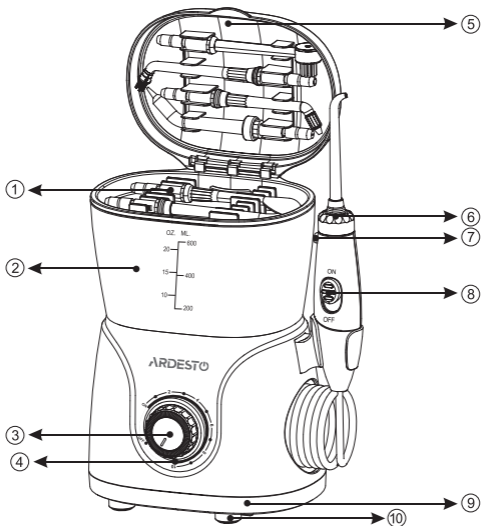PRODUCT STRUCTURE

1. Tips and storage box
2. Removable reservoir
3. OFF / ON
4. Water pressure regulator
5. Reservoir Lid box

6. 360° Rotating Tip
7. Tip Removing Knob
8. Water flow stop button
9. Base suction cups
10. Rubberized feet

PRODUCT DESCRIPTION

					
Ordinary High-Pressure Tip	Periodontal Tip	Orthodontic Tip	Dental Plaque Tip	Tongue Cleaner	Hydraulic rotating toothbrush

- Ordinary High-Pressure Tip: 3 color ring. For hygienic purpose. Recommend using by 3 people respectively (family pack).

The other 4 special tips

- Periodontal Tip x 1
- Orthodontic Tip x 1
- Dental Plaque Tip x 1
- Tongue Cleaner x 1
- Hydraulic rotating toothbrush x 1

See more instructions on page 22-23

This appliance should be kept out of reach of children at all times


Instructions

Plug into socket. If the socket is controlled by a wall switch, make sure the switch is turned on.

RECOMMENDED TECHNIQUE

Direct the jet stream at 90-degree angle to your gumline, slightly close lips to avoid splashing but allow water to flow freely from mouth into the sink. For best results, start in the molar area (back teeth) working towards the front teeth. Glide jet tip along gumline and pause briefly between teeth. Continue the process until all areas around and between the teeth have been cleaned.

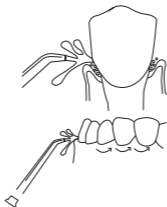
Pause Control

You can temporarily stop the flow at any time by slide down button on the handle.

When finished

Turn the unit off, any liquids left in the reservoir.
Unplug the unit.

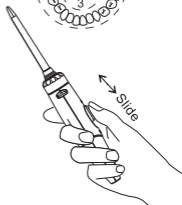Introduction to the jet tips

Periodontal Tip

The periodontal tip can direct water or liquid oral medicines in specific locations beneath the gum lining.

Ordinary High-Pressure Tip

Flosses deep locations within the oral cavity with a 360 degree range, supplementally covering places which are hard to clean, food scraps and harmful bacteria located deep within dental cracks or beneath the gum. This provides you with a pleasant dental experience.

Tongue Cleaner Tip

Place the tongue scraper on the tongue before switching on the Irrigator. Lightly place the tongue scraper on the centre-rear part of the tongue, lightly pulling it forwards while keeping it in contact with the surface of the tongue.

Dental Plaque Tip


Dental plaque is a group of bacteria that stick to the surface of the teeth yet has not mineralized. The jet tip is very effective at preventing dental plaque. The special-use dental plaque tip is for cleaning areas which is difficult to reach in braces and other dental devices.

Orthodontic Tip


The special Orthodontic Tip is uniquely designed to simultaneously brush and rinse hard-to-reach areas around braces and other dental work.

To use the Orthodontic Tip, lean low over the sink and use the brush bristles to clean around dental work. Water will flow through the brush to rinse away debris that is removed.

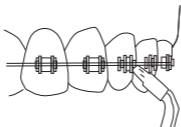
Hydraulic rotating toothbrush

The brush head is designed according to the dental scaling tools, which is different from the traditional toothbrush and just covers the whole tooth, thoroughly cleaning the surface of each tooth, the plaque hidden inside the tooth and the gingiva crevice, and will not hurt And gums.

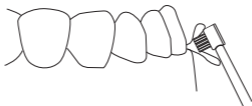
Preparing for reservoir

Remove the lid & reservoir from the base, Fill the reservoir with lukewarm water, Put the reservoir back to the base and press down firmly and put on the lid.

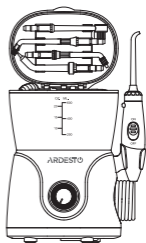
Inserting and removing tips

Insert the tip into the center of the knob at top of the dental the colored ring will be flush with the end of the knob if the correctly locked into place. To remove the tip from the handle, press the jet tip eject button and pull the tip from the handle.


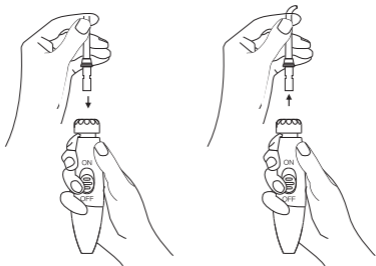
Adjusting The Pressure Setting

When you hear «click», said the switch has been opened. Then turn the clockwise to adjust the pressure from small to large. The machine will stop automatically in excess of 3 minutes at any stall. Next time when opening need to return to a «click», again to start the machine.

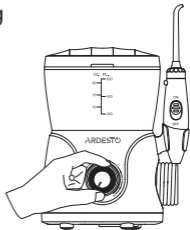Using the Irrigator

lean low over sink and place the jet tip in mouth. Aim the tip toward teeth. With your free hand, turn the on/off button on. Set pressure at desired setting. Direct the stream by turning the knob at the top of the handle.

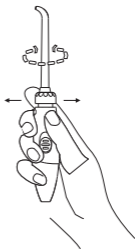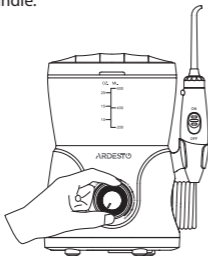CLEANING

After using any special solution in the reservoir, take a moment after every use to rinse the unit.

Partially fill the reservoir with warm tap water, aim the jet tip into your sink, turn on the unit and let it run until the reservoir is empty. This will help keep the tip from clogging and help your unit to last longer.

Clean the product when necessary by using a soft cloth and mild non-abrasive cleanser.

Before cleaning, please unplug from the electrical outlet.

Before exposing the dental jet to freezing temperatures, remove the reservoir and run the unit until completely empty.

SERVICE MAINTENANCE

1. After using the our irrigator, be sure to unplug the power cord, and the tank upside down on a flat surface, when the water is completely dry, put it into a cool, dry place together with the cleaned and dry host, it should be placed in the box or similar storage when not in use for a long time.
2. When a machine malfunctions, you can contact the seller to inquiry related maintenance methods.
3. Do not disassemble or attempt to resolve this product resulting in any problem, not within the company's warranty.

COMMON PROBLEMS AND SOLUTIONS

Problem	Cause	Solution
Jet Tips and Tongue Cleaner Leakage Between the top and handle knob	Tip not fully engaged	Remove tip and reinstall
Periodontal Tip Tear in soft tip	Tip is worn	Replace with new tip
Irrigator Base hose to the handle leaks	Damage to the hose	Can be replaced by the seller for a self-service repair kit
Reservoir leaks	Missing water tank valve or water tank valve is upside down	Replace Reservoir valve (contact the seller); Turn the valve right side up
Inadequate pressure	Reservoir not positioned properly	Be sure outlet is working, test with another appliance (check electrical reset button and or wall switch)
Unit without water flow out	Motor not running	Fill out reservoir

We warrant to the original purchaser/owner of this new product that it is free for defects in materials and workmanship for 1 year from the date of purchase.

CORRECT DISPOSAL OF THIS PRODUCT

This marking indicates that this product should not be disposed with other household wastes. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

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Warranty conditions:

The plant guarantees normal operation of the product within 12 months from the date of the sale, provided that the consumer complies with the rules of operation and care provided by this Manual. Service life of the product is 5 years.

In order to avoid any misunderstanding we kindly ask you to carefully study the Owner's Manual, the terms of warranty, and check the correctness of filling the Warranty Card. The Warranty Card is valid only if the following information is correctly and clearly stated: model, serial number of the product, date of sale, clear seals of the company-seller, signature of the buyer. The model and serial number of the product must correspond to those specified in the Warranty Card.

If these terms are violated, and if the information specified in the Warranty Card is changed, erased or rewritten, the Warranty Card is considered invalid. Setup and installation (assembly, connection, etc.) of the appliance are described in the attached documentation; they can be carried out both by the user, and by specialists having the necessary qualification sent by USC or selling companies (on a paid basis). In this case, the person (the institution) installing the appliance shall be responsible for the correctness and quality of installation (setup). Please pay attention to the importance of the correct installation of the appliance, both for its reliable operation, and for obtaining warranty and free service. Require the installation specialist to enter all the necessary information about the installation in the Warranty Card.

In case of failure of the appliance during the warranty period by the fault of the manufacturer, the owner is entitled to free warranty repair upon presentation of the correctly filled Warranty Card together with the appliance in the factory configuration packed in the original packaging to the warranty workshop or the place of its purchase. Satisfaction of the customer's claim due to the fault of the manufacturer shall be made in accordance with the law «On Protection of Consumers' Rights». In the case of warranty repair the warranty period shall be extended for the period of repair and shipment.

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WARRANTY CARD

ARDESTO™

Dear Buyer! Congratulations on your purchase of the ARDESTO brand appliance, which was designed and manufactured in accordance with the highest quality standards, and we thank you for choosing this particular instrument.

We ask you to keep the coupon during the warranty period. When purchasing a product, require a full warranty card.

1. Warranty service is carried out only if there is a correctly and clearly filled original warranty card, which indicates: product model, date of sale, serial number, warranty service period, and the seller's seal. *
2. The service life of household appliances is 5 years.
3. The product is intended for consumer use. When using the product in commercial activities, the seller / manufacturer does not bear warranty obligations, after-sales service is performed on a paid basis.
4. Warranty repair is carried out within the period specified in the warranty card for the product in an authorized service center on the conditions and terms determined by applicable law.
5. The product is withdrawn from the warranty in case of violation by the consumer of the operating rules set forth in the instruction manual.
6. The product is removed from warranty service in the following cases:
 - misuse and non-consumer use;
 - mechanical damage;
 - damage caused by the ingress of foreign objects, substances, liquids, insects;
 - damage caused by natural disasters (rain, wind, lightning, etc.), fire, domestic factors (excessive humidity, dust, aggressive environment, etc.)
 - damage caused by non-compliance of power and cable network parameters with state standards and other similar factors;
 - when operating equipment in the power supply network with a missing single ground loop;
 - in case of violation of seals installed on the product;
 - lack of serial number of the device, or inability to identify it.
7. The warranty does not cover consumables and accessories, as well as filters, shelves, drawers, stands and other containers for storing products.
8. The warranty period is 12 months from the date of sale.

* Tear-off maintenance tickets are provided by an authorized service center.

The completeness of the product is checked. I have read the terms of the warranty service, no complaints.

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