**ARDESTO™**

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## HEALTHY TEETH AND GUMS! EASY CLEANING OF ORTHODONTIC CONSTRUCTIONS!


You can effectively and easily take care of your gums, teeth, and beauty of your smile, using Ardesto irrigator on a daily basis! With the combination of water pressure and pulsation, the irrigator removes harmful bacteria and food debris deep between teeth and below gumline braces and other orthodontic constructions where brushing and even dental floss won't reach.

### Traditional Clean

Traditional brushing can't easily reach hidden food debris and harmful bacteria



Flossing has limited access below the gumline and may damage gums
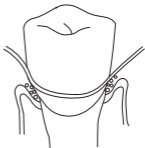
### Irrigator Clean

Removes trapped food debris and harmful bacteria with high frequency pulse water floss.

If you use irrigator on a daily basis, it gives you fresh breath, makes gums and teeth healthier, and helps to clean orthodontic constructions! Consult with your dentist on the irrigator modes usage.
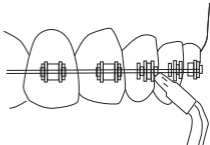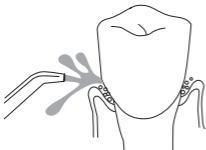
## CONTENTS

---

Important safeguards .....	21
Equipment .....	25
Product Accessories .....	26
Product Description .....	27
Operation Steps.....	28
Using Tips .....	28
Disassembly & Assembly .....	29
Cleaning .....	30
Service Maintenance.....	31
Battery Care and Disposal .....	34
Common issues .....	34
Correct disposal of this product.....	35
Warranty Conditions .....	38

## **IMPORTANT SAFEGUARDS**

---

When using electrical products, especially when children are present, the basic safety precautions below should always be followed.

### **READ ALL INSTRUCTIONS BEFORE USING.**

#### **DANGERS:**

- To reduce the risk of electrocution:
- Do not handle charger with wet hands.
- Do not immerse in water or other liquid.
- Do not use while bathing.
- Do not place or store product where it can fall or be pulled into a tub, shower stall or sink.
- Do not reach for a product that has fallen into water. Unplug immediately.

- Check the charger cord for damage before the first use and during the life of the product.

**To reduce the risk of burns, electric shock, fire, serious injury or damage issues:**

Do not plug this device into a voltage system that is different from the voltage system specified on the device or charger.

Do not use this product if it has a broken cord or plug. In case of not working properly, or it has been dropped into any liquid, contact your local distributor.

Do not direct water flow under the tongue, into the ear, nose or other delicate organs. This product is capable of producing pressures that may cause serious damage to these organs. Refer to operating instructions for correct usage.

Use this product only as

indicated in these instructions or suggestions by your dental professional.

It's normal phenomenon if gum bleed a little when using water flosser for the first time. If bleed much, please stop using and ask for your dentist advice.

Fill reservoir with water or other solutions recommended by dental professional.

Do not drop or insert any foreign object into the product.

Keep charger and cord away from heated surfaces.

Do not operate where oxygen or aerosol sprays are being used.

Do not use mouth washes/rinses that contain iodine, bleach, or tea treeoil.

Do not use while wearing any oral jewelry. Please put off prior to use. Do not use if you have an open

wound on your tongue or in your mouth.

If your physician advise you to receive antibiotic premedication, you should consult your dentist before using this instrument or any other oral hygiene aid.

For household use only. Do not use outdoors.

Closely instruct and supervise children and individuals with special needs in the proper use of this product.

Children should be supervised to ensure that they do not play with the appliance.

This device can be used by children over 8 years of age and persons with physical or sensory or mental disabilities or with insufficient experience and knowledge, if they are

under constant supervision or instructed in the safe use of the device and they understand the possible dangers. Cleaning and maintenance by the consumer should not be carried out by children without supervision.

## **EQUIPMENT**

---

Main body – 1 pc.

Water tank – 1 pc.


Jet tips – 2 pcs.

USB cable – 1 pc.

User manual with warranty card – 1 pc.

## PRODUCT ACCESSORIES

---

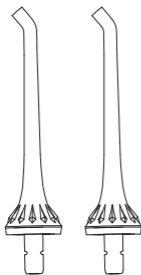Tips x 2


USB cable x 1

Replaceable tip/accessories are available from retail stores

## PRODUCT DESCRIPTION

---

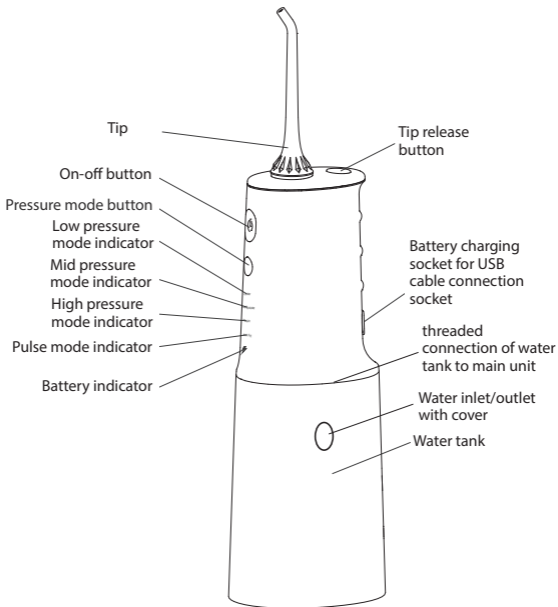
## OPERATING STEPS

### Start

**IMPORTANT!** Please fully charge the Water flosser for the first use by plugging the USB cable into the charging socket with power. The product cannot be used while charging.

USB adapter is not included. The USB cable is compatible with any 220V/5V, 1A USB power adapter.

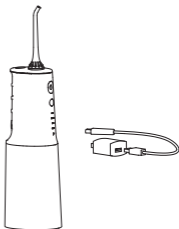
### USING TIPS

Tilt down to place the tip in the mouth, and point the tip at the tooth, then turn on the device.

Make water jet at 90° to the gingiva, and close your lips slightly to avoid water spray. For best results, start with the molars (back-teeth) and clean towards the front teeth. Slide the tip along the gingiva line and rest briefly between teeth until it is cleaned around gingiva and between teeth.
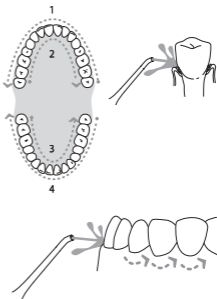
## DISASSEMBLY & ASSEMBLY


### After-use

#### Drainage of residual liquid

Turn off the device, and open the back tank lid with the injection hole facing down (shown in Fig. 1), then turn on the device after clearing in the internal channels of the irrigator the residual liquid in the tank, to drain the residual liquid in the until there is no more water to drain, then turn off the device.

Separate water tank from the main unit of the irrigator to wash or dry it. It's easy to do. Hold the product body with your left hand and hold the tank with the right hand, as shown in Fig. 2, and turn it counterclockwise to pull down the water tank, as shown in Fig. 3.
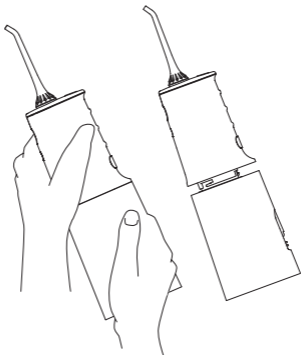
## CLEANING

---

After using any special solution in the reservoir, take a moment after every use to rinse the unit. Partially fill the reservoir with warm tap water, aim the jet tip into your sink, turn on the unit and let it run until the reservoir is empty. This will help keep the tip from clogging and help your unit to last longer.

Clean the product when necessary by using a soft cloth and mild non-abrasive cleanser. Before cleaning, please unplug from the electrical outlet.

Before exposing the irrigator to freezing temperatures, remove the reservoir and run the unit until completely empty.

## SERVICE MAINTENANCE

---


1. After using the irrigator, be sure to unplug the power cord, and the tank upside down on a flat surface, when the water is completely dry, put it into a cool, dry place together with the cleaned and dry host, it should be placed in the box or similar storage when not in use for a long time.
2. When a machine malfunctions, you can contact the seller to inquiry related maintenance methods.
3. Do not disassemble or attempt to resolve this product resulting in any problem, not within the company's warranty.

## Installation and removal of tips


Insert the nozzle into the middle hole of the tip knob until it clicks.

As shown in Fig. 4, open the tank lid and fill water into the tank; Note: Make sure that the device is off while filling the tank.

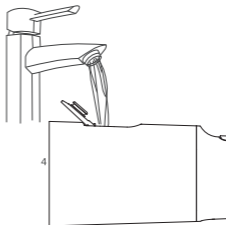If the tip needs replacing or removing, press the release button and while remove the tip.

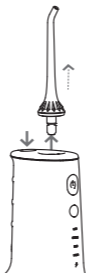- This product has memory function, so default work as the mode of last using.
- Press the mode button in the power-on state to switch between the low, mid, high and pulse modes; the working gear position and the power display are always on, and the other non-working gears are off.
- Mode selection can also be made when power off, and you can choose your favorite mode and then turn it on. Running in any mode for more than 120 seconds, the product will shut down automatically.
- When the product is charging, the power indicator light flashes; when the power indicator light is on, the product is fully charged.
- During use, when the power indicator light changes from white to red, the product is low battery; when the red light is flashing, the product needs charging.
- This product cannot be used while charging.
- Rotate the tip knob to adjust the direction of the water flow.
- Because gingiva bearing capacity differs from people to people, the first use of the Water Flosser may be accompanied with bleeding gingiva. This is normal. Please feel free to use it!
- In the first use of the portable water flosser, please start from the lowest mode and gradually shift to the higher with oral adaptation.

## BATTERY CARE AND DISPOSAL

---

### Battery Maintenance

- Lithium-ion battery: 3.7V, 1900mAh One 4-hour charge will be enough for a week. This type of battery is allowed for air transportation in your checked-in or cabin baggage. Please contact your air carrier for the transportation conditions.
- The device must be fully charged for 24 hours for the first use, since then charging about 4 hours.
- Plug the USB cable into the charging socket and power on;
- If it is idle for a long time, fully charge the device;
- Charge immediately once low-battery warning so that to utilize battery well.

### Disposing of the Battery

The device is equipped with a lithium battery, which will be recycled at the end of life.

Please take batteries, which are non-replaceable, to your local recycling center.

### FREQUENTLY ASKED QUESTIONS

---

Issues	Reason	Treatment
Leakage between the tip and the station	Broke seal ring in the tip holder and no tip in place	Contact the seller to repair the tip into place
Water tank leakage	Damaged seal ring	Contact the seller to repair or replace the seal ring
Water tank leakage	Low battery mode	Full charge
Failure to start	Low battery mode	Full charge
Full charge	Low battery mode	Cleaning of filter screen

If the above solutions do not work, please contact the sales store or the designated after-sales service center.

## CORRECT DISPOSAL OF THIS PRODUCT

---

This marking indicates that this product should not be disposed with other household wastes. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources.

To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

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## Warranty conditions:

---

The plant guarantees normal operation of the product within 12 months from the date of the sale, provided that the consumer complies with the rules of operation and care provided by this Manual. Service life of the product is 5 years.

In order to avoid any misunderstanding we kindly ask you to carefully study the Owner's Manual, the terms of warranty, and check the correctness of filling the Warranty Card. The Warranty Card is valid only if the following information is correctly and clearly stated: model, serial number of the product, date of sale, clear seals of the company-seller, signature of the buyer. The model and serial number of the product must correspond to those specified in the Warranty Card.

If these terms are violated, and if the information specified in the Warranty Card is changed, erased or rewritten, the Warranty Card is considered invalid. Setup and installation (assembly, connection, etc.) of the appliance are described in the attached documentation; they can be carried out both by the user, and by specialists having the necessary qualification sent by USC or selling companies (on a paid basis). In this case, the person (the institution) installing the appliance shall be responsible for the correctness and quality of installation (setup). Please pay attention to the importance of the correct installation of the appliance, both for its reliable operation, and for obtaining warranty and free service. Require the installation specialist to enter all the necessary information about the installation in the Warranty Card.

In case of failure of the appliance during the warranty period by the fault of the manufacturer, the owner is entitled to free warranty repair upon presentation of the correctly filled Warranty Card together with the appliance in the factory configuration packed in the original packaging to the warranty workshop or the place of its purchase. Satisfaction of the customer's claim due to the fault of the manufacturer shall be made in accordance with the law «On Protection of Consumers' Rights». In the case of warranty repair the warranty period shall be extended for the period of repair and shipment.

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# WARRANTY CARD

**ARDESTO™**

Dear Buyer! Congratulations on your purchase of the ARDESTO brand appliance, which was designed and manufactured in accordance with the highest quality standards, and we thank you for choosing this particular instrument.

**We ask you to keep the coupon during the warranty period. When purchasing a product, require a full warranty card.**

1. Warranty service is carried out only if there is a correctly and clearly filled original warranty card, which indicates: product model, date of sale, serial number, warranty service period, and the seller's seal. \*
2. The service life of household appliances is 5 years.
3. The product is intended for consumer use. When using the product in commercial activities, the seller / manufacturer does not bear warranty obligations, after-sales service is performed on a paid basis.
4. Warranty repair is carried out within the period specified in the warranty card for the product in an authorized service center on the conditions and terms determined by applicable law.
5. The product is withdrawn from the warranty in case of violation by the consumer of the operating rules set forth in the instruction manual.
6. The product is removed from warranty service in the following cases:
  - misuse and non-consumer use;
  - mechanical damage;
  - damage caused by the ingress of foreign objects, substances, liquids, insects;
  - damage caused by natural disasters (rain, wind, lightning, etc.), fire, domestic factors (excessive humidity, dust, aggressive environment, etc.)
  - damage caused by non-compliance of power and cable network parameters with state standards and other similar factors;
  - when operating equipment in the power supply network with a missing single ground loop;
  - in case of violation of seals installed on the product;
  - lack of serial number of the device, or inability to identify it.
7. The warranty does not cover consumables and accessories, as well as filters, shelves, drawers, stands and other containers for storing products.
8. The warranty period is 12 months from the date of sale.

\* Tear-off maintenance tickets are provided by an authorized service center.

**The completeness of the product is checked. I have read the terms of the warranty service, no complaints.**

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