**ARDESTO™**

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Before using the device – irrigator - please read this manual carefully and save it for future reference.

## **IMPORTANT SAFETY PRECAUTIONS**

- Use the device according to the instructions below or as recommended by your dentist.
- This appliance can be used by children from 8 years and above and people with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- The device is intended for domestic use only as intended in accordance with this manual. Do not use the device outdoors. Do not use the device for commercial purposes.
- Do not use the device outdoors and in places with high oxygen content or where aerosols (sprays) are used.
- Before connecting the device to the electrical network, make sure that the voltage indicated in the technical characteristics on the label of the device corresponds to the voltage parameters of your electrical network.
- The device is not designed to be controlled by an external timer or a separate remote-control system.

- Do not use the device with a damaged power cord or plug, or if the appliance is not working properly or has been damaged. Contact an authorized service center for diagnostics or repair of the device.
- Check the power cord and accessories for damage before each use. Do not use them if they are damaged.
- If the power cord is damaged, it must be replaced by the manufacturer, service center or similarly qualified personnel in order to avoid a hazard.
- Do not connect the device to the mains or remove the plug from the outlet with wet hands.
- Avoid contact of the device and power cable with heat sources.
- Do not place the device under direct sunlight.
- Install the device on a flat, stable surface.
- Do not close the ventilation openings of the device.
- Do not place or store the device where it could fall or slide into a bathtub, shower or sink.
- Do not use the device while taking a bath.
- Do not immerse the power cord, plug or device body in water or any other liquid to avoid electric shock or damage to the device.
- If the device has fallen into water, unplug it immediately. Contact an authorized service center.
- Do not leave the device connected to the mains unattended.
- Use only accessories recommended by the manufacturer.

- Do not insert foreign objects into the device or accessories.
- Always unplug the device after use or before cleaning.
- Before filling the water tank disconnect the device from the electrical outlet. Do not exceed the maximum fill level.
- Fill the water tank only with purified, distilled water or a solution recommended by the dentist.
- Do not direct the water jet under the tongue, ears, nose or other sensitive areas. This device creates pressure that can cause serious damage to these areas.
- Do not use the device if there are open wounds on the tongue or in the mouth.
- If you need to be premedicated with antibiotics before dental procedures consult your dentist before using this device or any other oral hygiene product.

**EQUIPMENT:**

Irrigator – 1 pc

Classic nozzle – 4 pcs

Dental plaque nozzle – 1 pc

Periodontal nozzle – 1 pc

Tongue cleaner nozzle – 1 pc


Instruction manual (warranty card is included) – 1 pc

## DEVICE DESCRIPTION

Thanks to a powerful pulsating jet of water, the irrigator helps to clean the teeth and remove harmful bacteria from hard-to-reach places between the teeth and gums.

- Simple and effective care for the health of gums and teeth
- Improvement of blood circulation and strengthening of gums, prevention of bleeding gums
- Care of braces
- Care of implants, crowns and other orthodontic constructions
- Removal of bad breath
- More effective and safer than dental floss

The irrigator in combination with the toothbrush provides comprehensive care for the teeth and oral cavity.
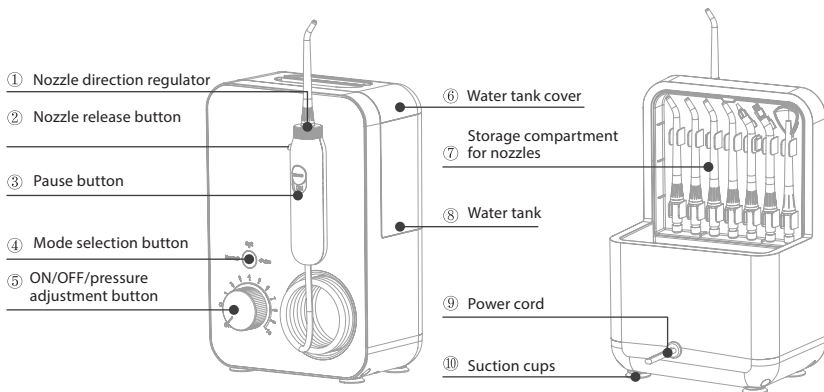
## **SPECIFICATIONS**

Removable water tank volume: 600 ml

Modes of water pressure regulation: 10

Pressure range: 30-125 PSI

Pulse frequency: 1400-1800 times per minute

Modes of care: normal, soft and pulse

3-minute smart timer

360° tip rotation

Electronics' waterproof class: IPX4

Voltage: 100-240 V~

Frequency: 50/60 Hz

Power: 20 W

**The manufacturer reserves the right to make changes to the technical specifications and design of products that do not impair performance without prior notice.**


## NOZZLE DESCRIPTION


Classic nozzle: completely cleans hard-to-reach areas in 360° range that toothbrush cannot reach.

Dental plaque nozzle: effectively prevents the formation of dental plaque and is intended for cleaning hard-to-reach areas, implants, dental bridges and crowns.

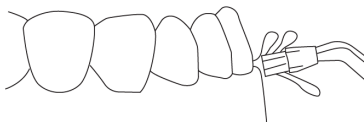Periodontal nozzle: massages gums, reduces periodontitis, gum disease and thickness of pockets.

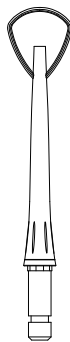
Tongue cleaner nozzle: effectively cleans tongue surface from harmful bacteria, keeps fresh breath.


## **DEVICE USAGE**

Remove and fill the water tank with water, then place it back into the device and press it lightly. Make sure the base of the water tank is aligned with the rubber seal.


Connect the device to the power source.

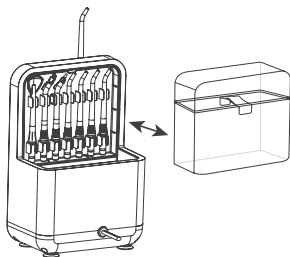Insert the nozzle into the handle until it clicks into place. To remove the nozzle press the nozzle release button and then remove the nozzle from the handle.

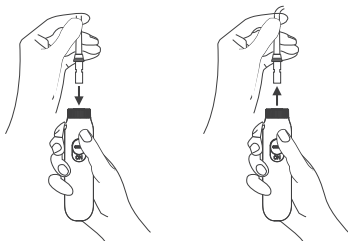
Place the spray head of the nozzle in the oral cavity and direct it to the teeth. Select the mode using the mode selection button. Start turning the switch to turn on the device and set the required pressure value. To increase the water pressure turn it clockwise.

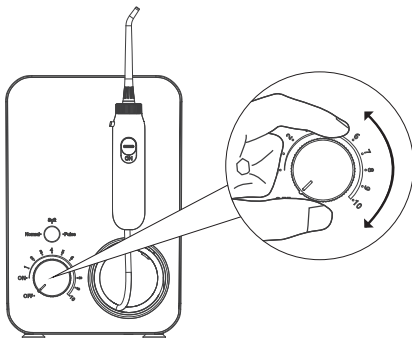
**Note:** gums should get used to the device, so the gums may bleed during the first uses. Use soft mode with minimal pressure first. In most cases the bleeding should stop in a week. If the condition persists for more than two weeks, stop using the device and consult your dentist.

The direction of the water jet can be adjusted by turning the nozzle direction regulator at the top of the handle.

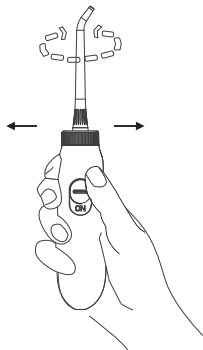**The device is equipped with a smart timer:** if the device is switched on for more than three minutes it will automatically turn off. To turn the device back on, turn the switch counterclockwise until the click is heard, then turn the switch clockwise.


**Pause button:** press the switch on the handle to stop water stream.
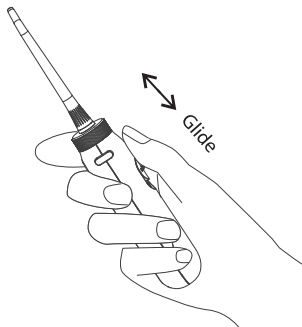
## Cleaning technique

### - Classic nozzle and dental plaque nozzle:

During irrigation the water jet should be directed at a right angle (90 degrees) to the surface of the teeth along the gum line. Cover the mouth to avoid splashing water, but do not cover it completely so that the water runs into the sink. For best cleaning results start brushing in the direction from the molars to the incisors. Move the spray head slowly along the gum line stopping periodically near interdental spaces.

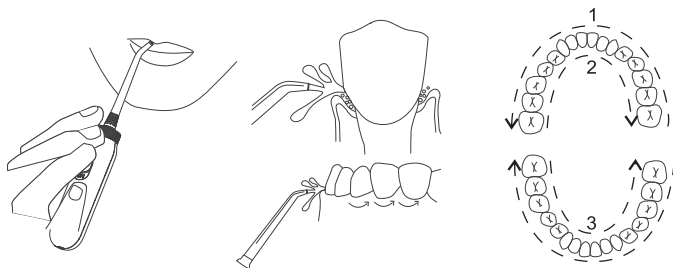**- Periodontal nozzle:**


Place the tip against a tooth at a 45-degree angle and gently place the tip under the gum line. Turn the device on to the lowest pressure setting and clean gently following gum line, stopping periodically near interdental spaces. Continue tracing along the gum line until finish of cleaning.

**- Tongue cleaner nozzle:**

Place the scraper in the center of the tongue, turn on the device and select a comfortable level of pressure. Gently move the scraper along the tongue without removing it from the surface of the tongue.

**After use** turn off the device, disconnect it from the mains and drain all water remaining in the tank.

**Storage compartment for nozzles:** for ease of use or cleaning the case with nozzles can be removed from the storage compartment by gently pulling out the buckle.
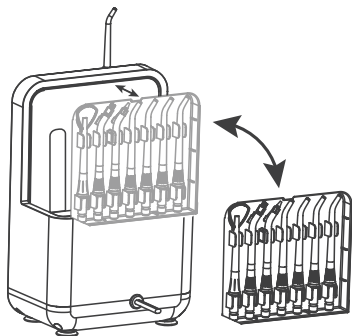
## **CLEANING AND CARE OF THE DEVICE**

- After each use and before cleaning switch off the device and remove the plug from the socket. Then remove the water tank and drain the water from it. Dry or allow the water tank to dry before installing back into the device.
- Press the nozzle release button and remove the nozzle, rinse it with water and then wipe it dry with a cloth. If the top part of the nozzle is deformed, replace the nozzle.
- Do not immerse the device in water or other liquids.
- Use a soft cloth and a neutral, non-abrasive detergent to clean the device. Wipe the outside of the device with a damp cloth, and then with a dry one.
- Do not use hard bristles, abrasive cleaners, solvents, etc. to clean the device, otherwise it may damage the device.
- Use a clean, dry cloth to clean the power cord.
- After using special solutions for the oral cavity rinse the device with clean water. Pour some warm water into the water tank and direct the spray head of the nozzle into the sink. Turn on the device and let it run until the water tank is empty. This will help avoid blocking of the spray head and extend the service life of the device.
- Do not store the device in places with very high or low temperature, direct sunlight and dusty environment. Store the device in a cool, dry place.

## **TROUBLESHOOTING**

<b>Malfunction</b>	<b>Possible cause</b>	<b>Solution</b>
Water leakage between nozzle and handle	The nozzle is installed incorrectly	Install the nozzle correctly
The tube that connects the device base with the handle leaks water	The tube is damaged	Contact an authorized service center
The water tank leaks	Water tank valve not installed or installed incorrectly	Install the valve correctly
Low water pressure	The water tank is installed incorrectly	Install the water tank correctly
The device does not turn on	Not properly connected to the power source or the engine does not work	Check whether the connection to the power source is correct or contact an authorized service center

**If the above actions did not help to eliminate the problem, do not try to repair the device yourself, contact an authorized service center.**

## **CORRECT DISPOSAL OF THIS PRODUCT**

Never dispose of the device with normal household waste. Dispose of the device through an authorized disposal company or an appropriate government agency. Please follow all applicable rules. If in doubt, contact waste disposal company.

Local community or municipal government can provide information on how to dispose of the product properly.

Dispose of the packaging materials environmentally. The packaging is made entirely of recyclable materials which should be disposed of at local recycling facilities.

Observe the marking of the packaging materials for waste separation, which are marked with abbreviations (a) and numbers (b) with following meaning:  
1–7: plastics / 20–22: paper and fiberboard / 80–98: composite materials.

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## Warranty conditions:

---

The plant guarantees normal operation of the product within 12 months from the date of the sale, provided that the consumer complies with the rules of operation and care provided by this Manual.

Service life of the product is 5 years.

In order to avoid any misunderstanding we kindly ask you to carefully study the Owner's Manual, the terms of warranty, and check the correctness of filling the Warranty Card. The Warranty Card is valid only if the following information is correctly and clearly stated: model, serial number of the product, date of sale, clear seals of the company-seller, signature of the buyer. The model and serial number of the product must correspond to those specified in the Warranty Card. If these terms are violated, and if the information specified in the Warranty Card is changed, erased or rewritten, the Warranty Card is considered invalid. Setup and installation (assembly, connection, etc.) of the appliance are described in the attached documentation; they can be carried out both by the user, and by specialists having the necessary qualification sent by USC or selling companies (on a paid basis). In this case, the person (the institution) installing the appliance shall be responsible for the correctness and quality of installation (setup). Please pay attention to the importance of the correct installation of the appliance, both for its reliable operation, and for obtaining warranty and free service. Require the installation specialist to enter all the necessary information about the installation in the Warranty Card.

In case of failure of the appliance during the warranty period by the fault of the manufacturer, the owner is entitled to free warranty repair upon presentation of the correctly filled Warranty Card together with the appliance in the factory configuration packed in the original packaging to the warranty workshop or the place of its purchase. Satisfaction of the customer's claim due to the fault of the manufacturer shall be made in accordance with the law «On Protection of Consumers' Rights». In the case of warranty repair the warranty period shall be extended for the period of repair and shipment.

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# WARRANTY CARD

Dear Buyer! Congratulations on your purchase of the ARDESTO brand appliance, which was designed and manufactured in accordance with the highest quality standards, and we thank you for choosing this product.

**Keep this card during the warranty period. When purchasing a product ask for the warranty card to be completely filled out.**

1. Warranty service is carried out only if there is a correctly and clearly filled original warranty card, which indicates: product model, date of sale, serial number, warranty service period, and the seller's seal. \*
2. The service life of household appliances is 5 years.
3. The product is intended for consumer use. When using the product in commercial activities, the seller / manufacturer does not bear warranty obligations, after-sales service is performed on a paid basis.
4. Warranty repair is carried out within the period specified in the warranty card for the product in an authorized service center on the conditions and terms determined by applicable law.
5. The product is withdrawn from the warranty in case of violation by the consumer of the operating rules set forth in the instruction manual.
6. The product is removed from warranty service in the following cases:
  - misuse and non-consumer use;
  - mechanical damage;
  - damage caused by the ingress of foreign objects, substances, liquids, insects;
  - damage caused by natural disasters (rain, wind, lightning, etc.), fire, domestic factors (excessive humidity, dust, aggressive environment, etc.)
  - damage caused by non-compliance of power and cable network parameters with state standards and other similar factors;
  - when operating equipment in the power supply network with a missing single ground loop;
  - in case of violation of seals installed on the product;
  - lack of serial number of the device, or inability to identify it.
7. The warranty does not cover consumables and accessories, as well as filters, shelves, drawers, stands and other containers for storing products.
8. The warranty period is 12 months from the date of sale.

\* Tear-off maintenance tickets are provided by an authorized service center.

**The completeness of the product is checked. I have read the terms of the warranty service, no complaints.**

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